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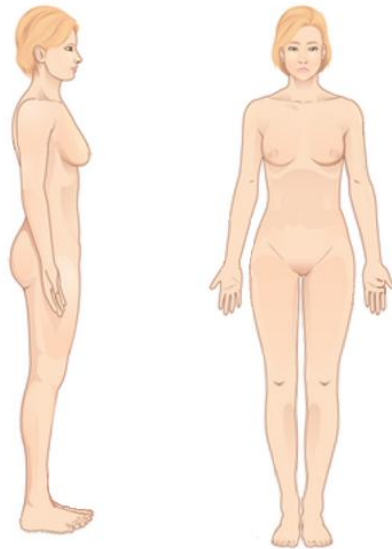
English :

Technical Requirements for Using the 3C Software

(3C Rapid Pain Relief 1.0)

With this software based on the “3C technique for rapid and precise pain relief,” you can immediately identify the point corresponding to the pain—that is, the point or area that needs to be massaged or needled.

When dealing with a patient in pain, accurate localization is essential before any treatment. To achieve precise localization, it is necessary to place the patient’s limbs in a standardized position.

Standardized position:

Both feet together, toes pointing forward, eyes looking straight ahead; both upper limbs relaxed along the sides of the body, palms facing forward, and the backs of the hands facing backward.

In this position, for the upper limbs:


- the radial side corresponds to the lateral (outer) side,
- the ulnar side corresponds to the medial (inner) side,
- the anterior surface is the flexor surface,
- the posterior surface is the extensor surface.

For the lower limbs:

- the fibular (peroneal) side corresponds to the lateral (outer) side,

- the tibial side corresponds to the medial (inner) side,
- the anterior surface is the extensor surface,
- the posterior surface is the flexor surface.

When identifying a painful point or searching for a corresponding point, you must strictly place the limb in this standardized position. For example, if the forearm is rotated (internal rotation), confusion will occur: at the wrist, the lateral side will shift toward the front or the medial side; at the elbow, internal rotation will also cause the lateral side to shift approximately halfway along the wrist pathway. In such cases, the localization of the painful point or corresponding point will inevitably be incorrect.

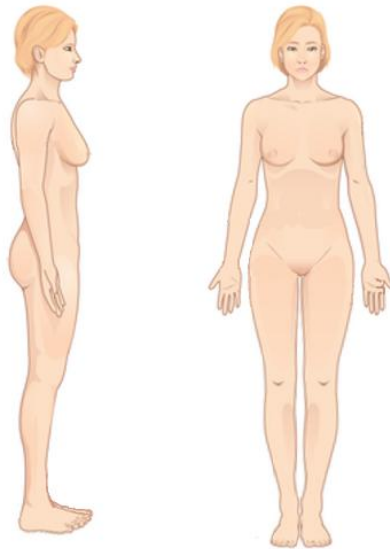After opening the software and selecting the language, you will see a measurement tool in the toolbar, specifically designed for the head and the thoracoabdominal region. After clicking the measurement tool, first click on point O, then move the mouse to the painful area; the red dot indicates the corresponding point.

For the limbs, you can directly click on the painful area; both the corresponding point and the painful area will appear in red.

The images do not indicate left or right, but arrows are shown at the bottom. The arrow points toward the side of the computer operator, and the image therefore represents the patient's side. If you need the opposite side, you can click the image flip button in the toolbar. This software was designed to enable as many people as possible to locate points accurately. As for the deeper theoretical foundations and operating mechanisms of the 3C technique, you may refer to our book *"Rapid Pain Relief with the 3C Technique."* This book is available in Chinese, French, and English. The purchase link can be found on our association's website:

<https://santemtc.com>

On our association's website, you will also find several introductory videos on the 3C technique to support your learning.

You may also join the "3C 2025 Sweden WeChat Learning Group." You can do so through a friend's recommendation, by sending us an email (contact@santemtc.com), by leaving a comment under the videos, or by writing to us under the section *"Testimonials from Acupuncturists after Learning the 3C Technique"*:

<https://santemtc.com/3c-temoin>

Thank you for your trust and support.

International Natural Health Association

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